

NEW!

Restoratives



Demi™ Pro Curing Light

The Light That Moves With You



All-Day Comfort

Reduce fatigue with the lightweight, ergonomic design.



Deeper Cures

Delivers consistent curing performance at clinical distances, thanks to its collimated output and improved beam profile.



Easy Access

Reaches every angle with the unique 360° and low profile rotating tip.



Safety

Minimize the risk of overheating with patented PLS (Periodic Level Shift) Technology and an integrated radiometer.

Cure With Control, Not Fatigue

- ✓ Ergonomic comfort
- ✓ Pen-like precision in a lightweight, balanced body at just 115g
- ✓ 360° access with a low-profile tip



360°
PRECISION
REACH EVERY ANGLE


"The Demi Pro curing light is very lightweight and made curing at every angle comfortable."



- Dr. Shahira Saad

Deeper Cures. Powerful Results.

Delivers powerful, collimated light output even at increased distances, ensuring reliable curing performance across various clinical scenarios.

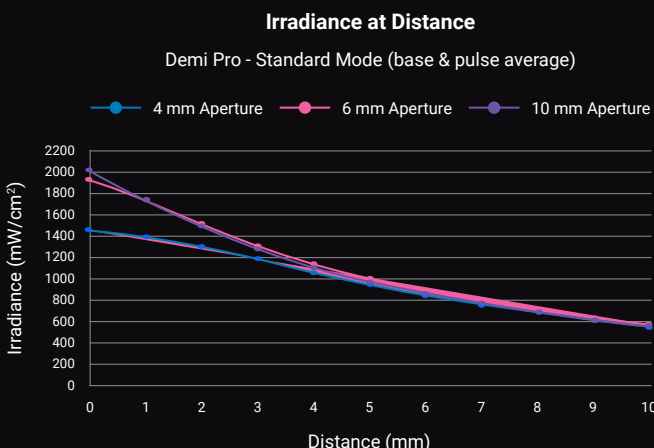

Deeper Cures With Consistency

Demi Pro delivers reliable curing results with consistent, high-intensity light distribution at clinical distances.


	Demi Pro	VALO™ X	SmartLite® Pro
Mode	Standard	Std Power (Steady Blue)	Standard
Test Result (0mm) Light tested directly on sensor	1393 mW/cm²	1137 mW/cm²	1233 mW/cm²
Test Result (6mm) Light tested using tooth	1025 mW/cm²	830 mW/cm²	748 mW/cm²

Bluelight Analytics Inc. - Demi Pro Evaluation Final Report | May 2026

Confidence Starts With Safety

Demi Pro includes a built-in radiometer to easily test its output, giving you peace of mind going into every procedure.

Demi Pro employs the power of Periodic Level Shifting, (PLS), shifting output intensity multiple times throughout the curing cycle to drive continuous high output without generating excess heat.

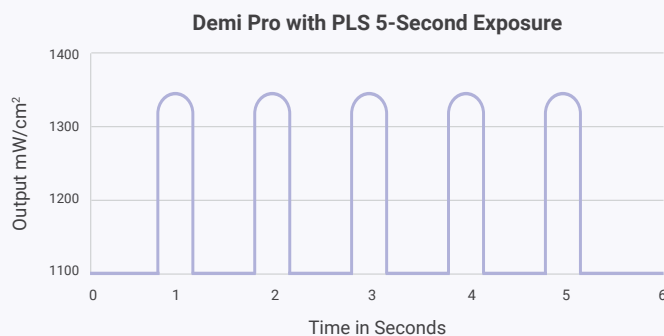

Transillumination Tip*

A powerful tool to aid in intraoral visualization, diagnosis, and treatment planning.



Ordering Information

Kit

37380 Complete Kit
1x Curing Tip, 1x Handpiece,
1x Charging Base, 2x Batteries,
3x Light Shields, 1x Tack Cure
Lens, 1x AC Adapter

Accessories

37381 Curing Tip
37382 Handpiece
37383 Charging Base
37384 Battery
37385 Light Shield Refill
37386 Tack Cure Lens
37387 Transillumination Tip
37379 AC Adapter

*Sold separately

